

DXL IMPLANT PLUS

CONICAL CONNECTION

Conical connection with inner hex provides precision implant-abutment connection and therefore, influences well on the stability process of permanent dental prostheses.

Conical connection is absolutely bacterial dense and shows excellent results in compression reduction in implant cervical part.

THREAD

Spiral-shaped implant with regular thread was develop for a better primary stability, therefore it is recommended to be used in D1, D2, D3 and D4 bone types.


Conical implant body perfectly complies with the drill shape for more precision bone fitting. Implant body tapering ensures better self-tapping property. Noninvasive apical implant part prevents from anatomical damaging during sinus lifting

SURFACE


DXL implants need a specific surface roughness to be incorporated fast, permanently and without complications into the jaw tissue. DXL Implants SLA surface is the latest generation of dental implant surfaces. The increased surface area provides perfect bone implant connection. Due to this advantage the implant has higher primary stability and quick osseointegration. It can be used with one- or two-stage protocol.

APICAL PART

Due to the precise design of implant apical part it can prevent from damage of anatomical bone structure and can give the opportunity to the implantologist to decrease the time of surgical implantation stage.

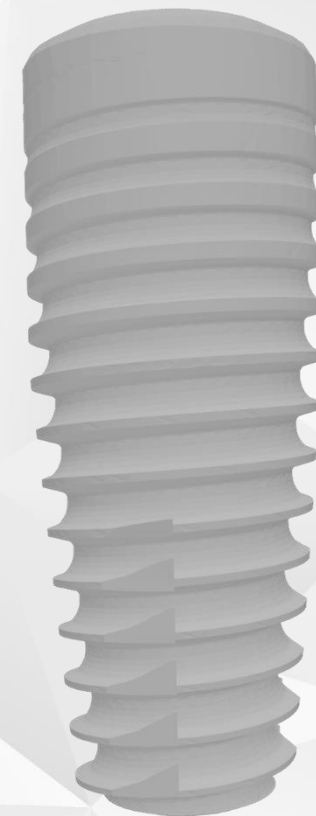

DXL IMPLANT PLUS

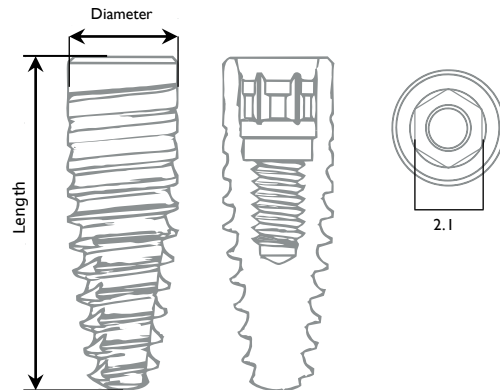

PRODUCT CODE
DIAMETER

DXL-A330080

8 mm



DXL-A330100

10mm



DXL-A330125

12,5mm

NARROW
PLATFORM

DIAMETER

3.3

COLOUR CODE

PINK

INTERNAL HEX

2.1

It is possible to order implants with cover screw Add -C to the relevant article and Add with cover to the relevant implant name in your order.

ATTENTION!

It is used for installation in front and lateral parts canines, incisors, premolars .

Please, consult DXL web site for Implantation Protocol.




www.dxlimplant.de / ifu.dxlimplant.de

DXL NARROW COVER SCREW STERILE
Ø2,6xH0mm / DXL-3000001-N-S



Only DXL Premium Hybrid
Implant with cover series
include cover screw.

DXL IMPLANT PLUS

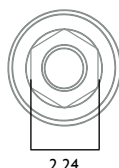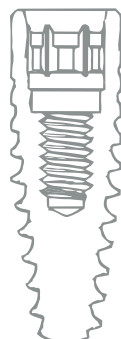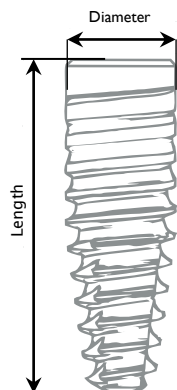

PRODUCT CODE
DIAMETER

DXL-A370080

8 mm

DXL-A370100

10 mm

DXL-A370125

12.5mm

DXL-A370150

15mm

STANDARD

PLATFORM

DIAMETER

3.75

COLOUR CODE

RED

INTERNAL HEX

2.24

It is possible to order implants with cover screw Add -C to the relevant article and Add with cover to the relevant implant name in your order.





DXLCOVERSCREWSTERILE
Ø3,1xH0mm/ DXL-3000001-S

Only DXL Premium Hybrid
Implant with cover series
include cover screw.

Please, consult DXL web site for Implantation Protocol.

www.dxlimplant.de / ifu.dxlimplant.de

DXL IMPLANT PLUS

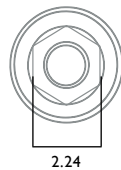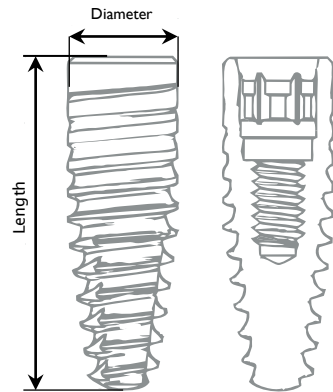

PRODUCT CODE
DIAMETER

DXL-PH-4108

8 mm

DXL-PH-4110

10 mm

DXL-PH-4112

12.5mm

DXL-PH-4114

15mm

STANDARD
PLATFORM

DIAMETER

4.1

COLOUR CODE

GREEN

INTERNAL HEX

2.24

It is possible to order implants with cover screw Add -C to the relevant article and Add with cover to the relevant implant name in your order.




DXL COVER SCREW STERILE
Ø3,1xH0mm / DXL-3000001-S

Only DXL Premium Hybrid
Implant with cover series include
cover screw.

Please, consult DXL web site for Implantation Protocol.

www.dxlimplant.de / ifu.dxlimplant.de

DXL IMPLANT PLUS

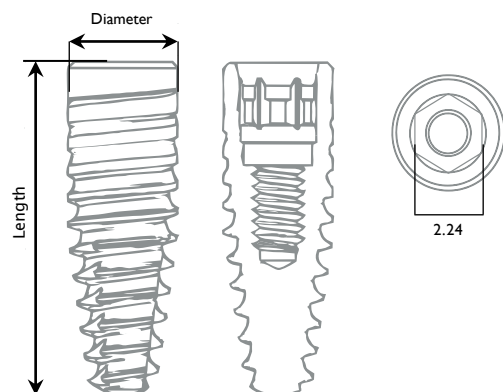

PRODUCT CODE
DIAMETER

DXL-PH-4I108

8 mm

DXL-PH-4I110

10 mm

DXL-PH-4I112

12.5mm

DXL-PH-4I114

15mm

STANDARD

PLATFORM

DIAMETER

4.8

COLOUR CODE

BLUE

INTERNAL HEX

2.24

It is possible to order implants with cover screw Add -C to the relevant article and Add with cover to the relevant implant name in your order.



DXL COVER SCREW STERILE
Ø3,1xH0mm / DXL-3000001-S

Only DXL Premium Hybrid
Implant with cover series include
cover screw.

Please, consult DXL web site for Implantation Protocol.

www.dxlimplant.de / ifu.dxlimplant.de